


# Materials Processing Data Sheet

## Microfluidizer® Processor Laboratory Results

Industry:	Chemical
Product:	10% suspension of Germanium, GeO <sub>2</sub> , in isopropyl alcohol
Test Objective:	Reduce germanium particle size from 4-5 µm to less than 1 µm
Processor:	M-110Y
Chambers:	<ul style="list-style-type: none"><li>• Diamond Interaction Chamber (upstream) G10Z (87 µm)</li><li>• Ceramic Auxiliary Processing Module (downstream) H30Z (200 µm)</li></ul>
Number of passes:	1
Shear per pass:	5,800,000 sec <sup>-1</sup>
Process Pressure:	23,000 psi
Particle Size:	Premix: 8.301 µm mean particle size [D <sub>90</sub> = 14.469 µm] 1 Pass: 0.502 µm mean particle size [D <sub>90</sub> = 0.612 µm]

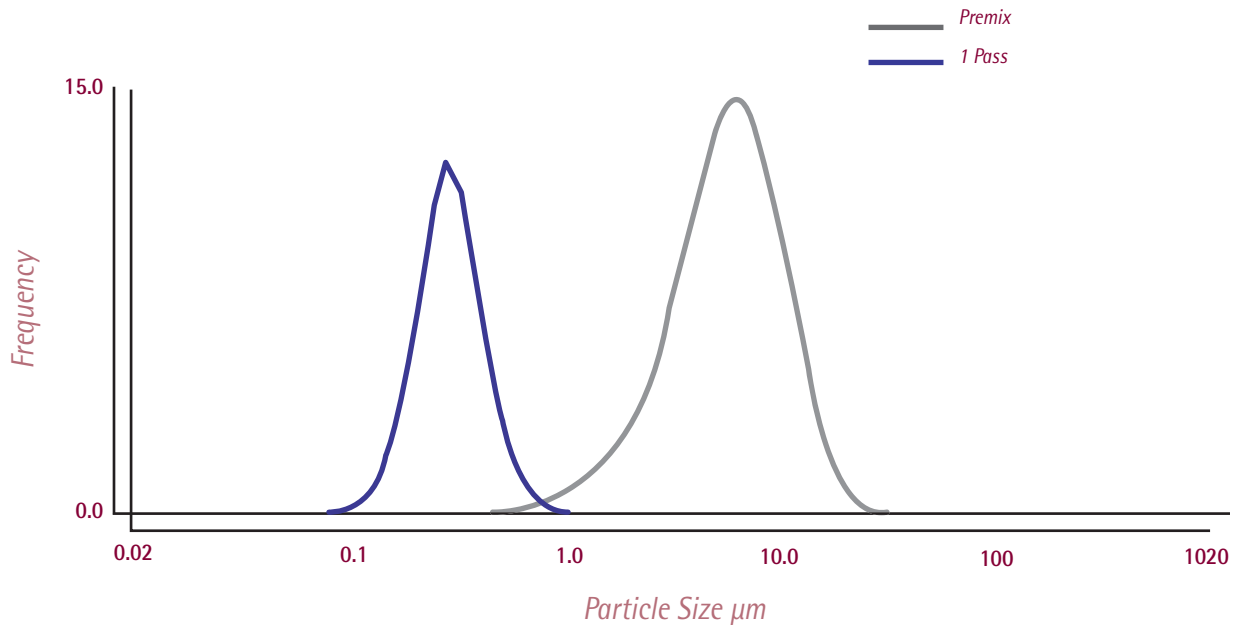

### Summary of Results:

The mean particle size of the germanium was reduced from 8.301 µm to 0.502 µm after one pass @ 23 kpsi using the M-110Y Microfluidizer processor.

# Materials Processing Data Sheet

## HC-2000 GMP Pneumatic

### HC-5000 GMP Pneumatic

The HC-2000 is our lightweight, easy to use, low pressure homogenizer. With sanitary fittings and 300 series stainless steel parts, processing with the HC-2000 meets GMP requirements. The HC-2000 is suitable for batch and continuous processing with flowrates to 2 liters/minute. The HC-5000 is identical to the HC-2000 but with higher pressure capabilities (to 5,000 psi) and flow rates to 950 ml/minute.

Pressure Range 250–2,000 psi

Pressure Range 500–5,000 psi

## HC-8000 Sanitary Pneumatic

The HC-8000 complies with sanitary standards to meet demanding product specifications and tight quality control. This unit employs a flow-through, sanitary pressure transducer with digital display and generates liquid pressures up to 8,000 psi and flowrates to 1 liter/minute.

Pressure Range 1,000–8,000 psi

## M-110S Small Volume Pneumatic

The M-110S is an air-powered unit designed to process smaller sample volumes than the other Microfluidizer models with greater product recovery. For samples as small as 14 ml (with >12ml recovery) up to continuous operation.

Pressure Range 3,000–23,000 psi

## M-110Y High Pressure Pneumatic

The M-110Y provides the highest operating pressure of our laboratory units. It is specifically designed for applications that require optimum pressure, such as cell rupture. The M-110Y is ideally suited for improving product yield and minimizing downtime.

Pressure Range 3,000–23,000 psi

## M-110EH Electric-Hydraulic

The M-110 EH is specifically designed for those customers that prefer an electrically-driven unit over a pneumatically-driven one. The unit reaches pressures as high as 25,000 psi and is mounted on locking casters for mobility.

Pressure Range 2,500 - 25,000 psi

## M-140K High Pressure Electric-Hydraulic

The M-140K is built to sustain prolonged high pressure demands and provide users with the opportunity to create new formulations, and improve existing products. The unit reaches pressures as high as 40,000 psi and is mounted on locking casters for mobility.

Pressure Range 8,000 - 40,000 psi

## M-210EH Electric-Hydraulic

The M-210 EH Series is the pilot plant line specifically designed to utilize an electric-hydraulic pump in place of the standard air-powered pump. Results achieved on laboratory Microfluidizers will scale up to this level. These units can achieve pressures up to 30,000 psi and flowrates to 1.12 gpm. The M-210EH is mounted on locking casters for mobility.

Pressure Range 2,500–30,000 psi

## M-700 Series

The M-700 Series offers pilot and production performance levels with specific features amenable to industrial applications, including coatings, waxes, inks, pigments, sealants, etc. The units are available in 15, 25, and 50 horsepower models, and can be fitted with 1 or 2 intensifier pumps that can be operated with independent product streams and flow rates.

Pressure Range to 40,000 psi

## M-610 Custom Designed Electric-Hydraulic Production

The M-610 Series Production Microfluidizer is capable of flowrates to 18 gpm and can process samples from one gallon to continuous operation. These custom built units employ an electric/hydraulic intensifier pump which can be designed to generate process stream pressures up to 40,000 psi.

Pressure Range to 40,000 psi



HC-8000



M-700



M-610

Microfluidics International Corporation

90 Glacier Drive, Suite 1000 • Westwood, MA 02090, USA

Tel: 617-969-5452 • 800-370-5452 • Fax: 617-965-1213

Email: [mixinginfo@idexcorp.com](mailto:mixinginfo@idexcorp.com) • [www.microfluidicscorp.com](http://www.microfluidicscorp.com)

IDEX Material Processing Technologies

Bramley Drive, Vale Park West, Evesham, Worcestershire, WR11 1JH, UK

Tel: (+44) (0) 1386 769 007 • Fax: (+44) (0) 870 1911116

Email: [mixinginfo@idexcorp.com](mailto:mixinginfo@idexcorp.com) • [www.microfluidicscorp.com](http://www.microfluidicscorp.com)



Microfluidics

IDEX  
MATERIAL PROCESSING